



USER MANUAL HB-P01

- Please carefully read all the cautions and instructions before using the device.
- Service fee may apply for any damage caused by misuse from not following the addressed cautions.

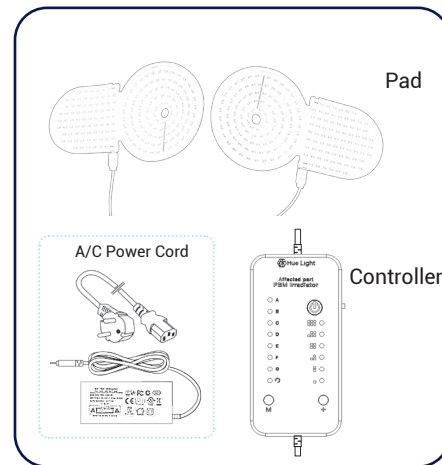


1 Start PBM Therapy

- 1 Connect the cords to connect the controller and the pad properly.
- 2 Connect the power cord to the proper outlet.
- 3 Place each pad inside the bra to cover the breast and the armpit.
- 4 Set the time using the switch on the side.
- 5 Set the light emission level by pressing the + button.
- 6 Set the Nogier Frequency Level by pressing the M button.
- 7 Press the Start button to start the session.

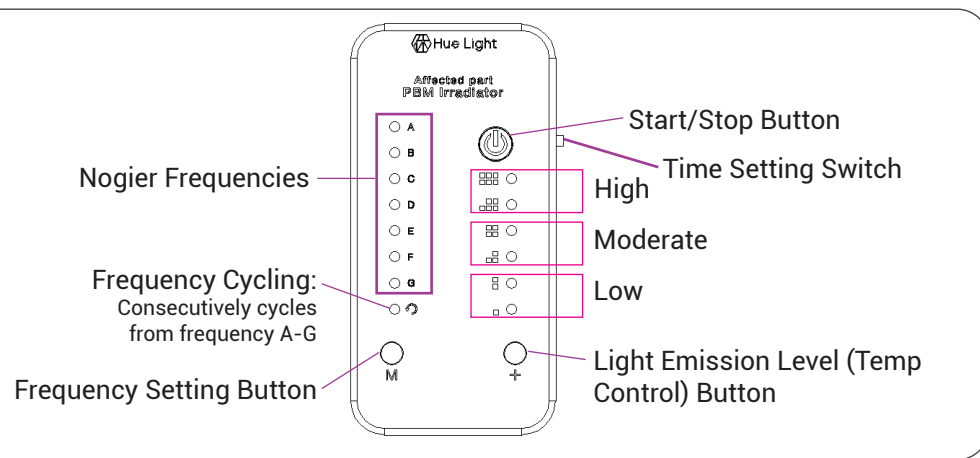Labels: Start/Stop Button, Time Setting Switch, High, Moderate, Low, Light Emission Level (Temp Control) Button, Frequency Cycling: Consecutively cycles from frequency A-G, Frequency Setting Button (M), Nogier Frequencies (A-G).



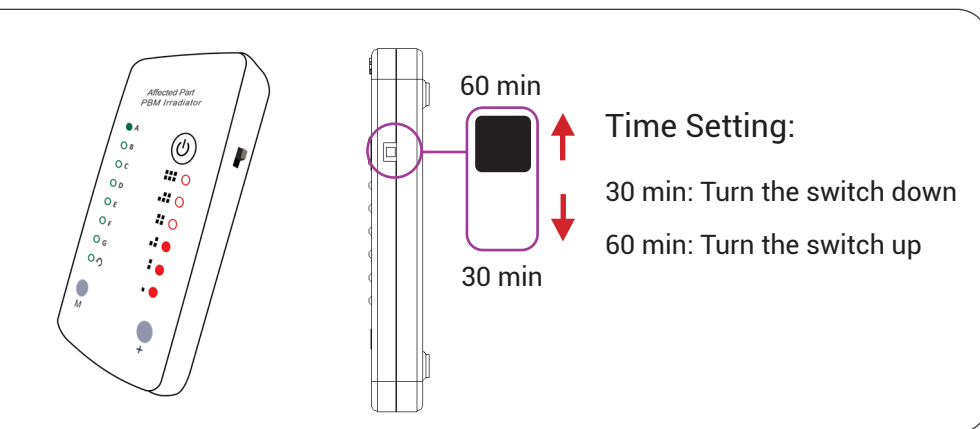

Time Setting:
30 min: Turn the switch down
60 min: Turn the switch up


! Cleaning & Storage

Cleaning: Use a tissue with alcohol or clean water to clean the pad after use. Wipe it dry with a dry cloth.
Storage: Disconnect all the cords and put the product inside the package box after use.

4 Nogier Frequencies



Frequency A	Wound healing Skin and nerve repair Reduce scar tissue, inflammation, and tumors
Frequency B	Nutritional assimilation Relieve allergy issues, edema, and lymph issues Balance nervous system
Frequency C	Muscle, skeletal, and myofascial pain relief
Frequency D	Balance the two sides of the brain Reduce and treat stress
Frequency E	Recover from spinal and skin disorders Pain control Reduce excess calcification (arthritis)
Frequency F	Recover from muscle spasms, facial pain, headaches, and depression Enhance circulation
Frequency G	Enhance memory Reduce psychological disorders and nervousness

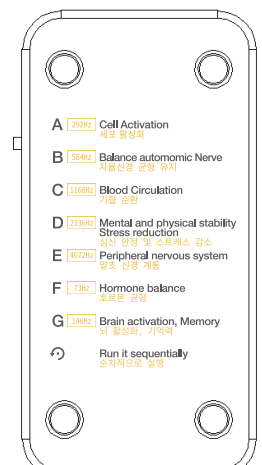

Safety and Maintenance

Cautions

- ◆ The highest heat level of the operating device is 38–40 °C. If you feel that the heat is too high, please lower the heat level by lowering the light emission level.
- ◆ The pad has LED chips inside the silicone. So, do not strike or put any object on the product. Piercing or folding the product too much may cause malfunctions.
- ◆ The product is not fully waterproof. Do not put it in water.

A/S Warranty

- ◆ If there is any issue with this product, you can get a free A/S service for 1 year. (A service fee may be applied for any cases that fall into the paid service category even during this warranty period.)
- ◆ **Paid Service**
 - Product failure or damage due to the user's misuse
 - Failure due to improper repair or modification
 - Failure caused by a natural disaster
 - Damage caused by failure to observe any cautions specified in this manual